

# GoodKnight™ H<sub>2</sub>O

HEATED HUMIDIFIER




- Effective humidification the whole night
- Easy to use
- Easy to clean

## Improve Comfort and Compliance for CPAP Users.

For CPAP patients experiencing upper airway dryness, Tyco Healthcare offers a mist of relief with the GoodKnight™ H<sub>2</sub>O Heated Humidifier.

Designed for use with the GoodKnight family, the GoodKnight H<sub>2</sub>O Heated Humidifier increases air humidity, reducing patient dryness in the nose and throat. By improving comfort, it helps promote patient compliance with sleep apnea treatment.

In addition, the humidifier features a small footprint that fits neatly with GoodKnight sleep therapy devices, creating the smallest, lightest CPAP/humidifier combination on the market.

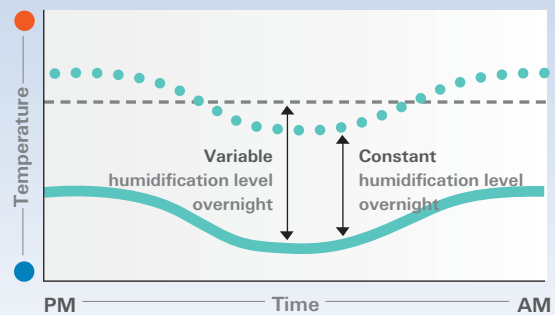

### HumidControl Feature

The HumidControl feature ensures a consistent temperature difference between the water and ambient air to minimize condensation in the tubing and keep the humidification level constant.

Room Temperature ———  
GoodKnight H<sub>2</sub>O .....  
Constant - - - - -  
Temperature Humidifier ———

**tyco**

Healthcare

**Puritan Bennett**

# GoodKnight H<sub>2</sub>O

HEATED HUMIDIFIER



GoodKnight™ H<sub>2</sub>O  
is compatible with:



GoodKnight™ 420 family

&



GoodKnight™ 425 Bi-Level  
devices

## Specifications

### Performance:

- Maximum flow: 120 liters per minute (leak flow included)
- Resistance: 0.2 cm H<sub>2</sub>O
- Service maximum pressure: 25 cm H<sub>2</sub>O
- Humidification rate: 12 mg H<sub>2</sub>O/L at 60 L/min and at the maximum setting
- Time of heating: ≤ 45 minutes

### Physical Characteristics

- Dimensions (W x H x L): 170 x 95 x 280 mm
- Weight (when empty): 0.8 kg
- Operational volume: 213 cm<sup>3</sup> (minimum volume of gas in the chamber)
- Volume usable of the water chamber: 350 ml
- Maximum level of liquid: indicated on the chamber

**Effective humidification:** The water chamber holds 350 ml of water that allows it to last the whole night.

**Easy to clean:** Comes with an easy-to-detach water chamber for trouble-free maintenance.

**Easy to use:** Features a rotary button that allows patients to adjust heating level.

**Space saving design:** Adds only 7.5cm length to the base of any GoodKnight 420/425 sleep therapy device.

## Ordering Information

Description	Catalog Number
GoodKnight H <sub>2</sub> O Heated Humidifier	M-114600-IN
GoodKnight H <sub>2</sub> O Water Chamber	M-314600-00
GoodKnight H <sub>2</sub> O Connection Adapter water chamber/CPAP	M-414600-05
GoodKnight 420G CPAP System	M-113900-IN
GoodKnight 420S CPAP System	M-113903-IN
GoodKnight 420E Auto-CPAP System	M-113902-IN
GoodKnight 425 Bi-Level Device	M-114500-IN
GoodKnight 425ST Bi-Level Device	M-114502-IN

Specifications subject to change without notice.

Nellcor Puritan Bennett Inc. © 2004 Tyco Healthcare Group LP. All rights reserved.

C-AD-GKH20/INT

**tyco**  
Healthcare

**Puritan Bennett**

Tyco Healthcare Deutschland GmbH  
Gewerbepark 1  
D-93333 NEUSTADT/DONAU  
Tel.: +49 9445 9590  
Fax: +49 9445 959155

Tyco Healthcare Norden AB  
Gårdsvägen 18  
S-16967 SOLNA  
Tel: +46 (0) 8585 60 505  
Fax: +46 (0) 8585 60 506

Tyco Healthcare Spain S.L.  
Avda. San Pablo, 28  
E-28820 COSLADA, MADRID  
Tel.: +34 91 669 68 48  
Fax: +34 91 674 03 24

Tyco Healthcare Italia SpA  
Via Rivoltana, 2/d  
I-20090 SEGRATE (Mi)  
Tel: +39 02 70317.1  
Fax: +39 02 7530125

Tyco Healthcare Nederland BV  
Hogeweg 105  
NL-5301 LL ZALTBOMMEL  
Tel: +31 (0) 418 57 66 00  
Fax: +31 (0) 418 57 67 93

Tyco Healthcare Belgium N.V.-S.A.  
Generaal de Wittelaan 9/5  
B-2800 MECHELEN  
Tel: +32 (0) 15 29 44 50  
Fax: +32 (0) 15 29 44 55

Tyco Healthcare Austria GmbH  
Campus 21, Europaring F09402  
A-2345 BRUNN/GEBIRGE  
Tel: +43 (0) 2236 37 88 39  
Fax: +43 (0) 2236 37 88 40

Tyco Healthcare Switzerland Ltd.  
Roosstrasse 53, Postfach  
CH-8832 WOLLERAU  
Tel: +41 1 786 50 50  
Fax: +41 1 786 50 10

Tyco Healthcare France S.A.  
3, avenue du Canada - Les Ulis  
F-91975 COURTABOEUF  
CEDEX  
Tel: +33 (0) 1 69 82 15 18  
Fax: +33 (0) 1 69 82 15 28

Tyco Healthcare Regional Export  
154 Fareham Road  
GOSPORT  
HAMPSHIRE  
UK - PO13 0AS  
Tel: +44 (0) 13 29 22 42 58  
Fax: +44 (0) 13 29 22 40 86

Tyco Healthcare Portugal  
Produtos de Saúde lda.  
Est: Outeiro Polima, Lote 10-1ºPiso  
Abóboda  
P-2785-521 S. Domingos de Rana  
Tel: +351 21 448 10 00  
Fax: +351 21 445 05 88

Tyco Healthcare South Africa  
Birchwood Court  
Montrose Road, P.O. Box 8108  
Halfway House  
1685 MIDRAND  
SOUTH AFRICA  
Tel: +27 (0)11 265 2000  
Fax: +27 (0)11 265 2050

Tyco Healthcare UK Limited  
154 Fareham Road  
GOSPORT  
HAMPSHIRE  
UK - PO13 0AS  
Tel: +44 (0) 13 29 22 41 59  
Fax: +44 (0) 13 29 22 40 83  
08/2004